

GELITA® EC

Enteric Capsule Release



Enteric Capsule Release GELITA® EC



- First of its kind
- Production with one gelatine product in one process
- Creates value-added market opportunities

GELITA® EC

Expanding GELITA's Portfolio of Capsule Release Possibilities

GELITA offers a comprehensive portfolio of gelatines that provides a broad release profile in soft capsules. Besides variety of versatile and tailor-made standard gelatines, there are GELITA® RXL and GELITA® RXL Advanced for improved capsule shelf-life and stability. These special gelatines provide reliable release and reduced cross-linking properties, even under extreme storage conditions. This GELITA capsule gelatine portfolio is complemented by GELITA® RXLR² for revolutionary fill release. The soft capsules produced with GELITA® RXLR² offer faster release of the fill and offer the added advantage of the reduced cross-linking properties of the RXL technology.


And, now – on the opposite end of the release profile spectrum, GELITA offers GELITA® EC; a gelatine product for the one-step production of truly enteric capsules.

Release Profile of the GELITA® Portfolio





Dissolution results

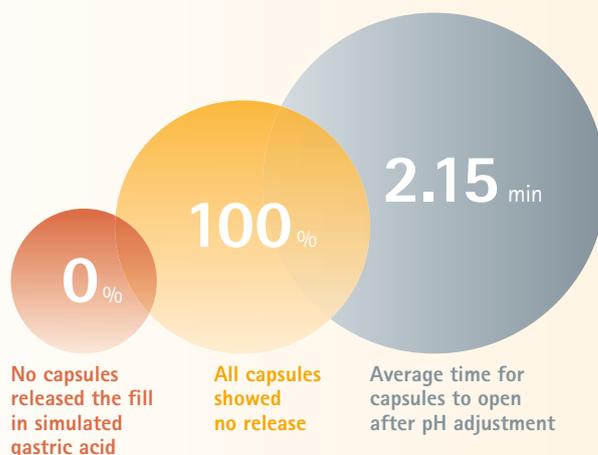

! Satisfying Consumer Demand

Consumers have been ingesting nutrients in soft capsules for decades due to their easy of swallowing and convenience. However, one nagging issue has been plaguing consumers for as long as they've been ingesting fish oil and other omega-3 fatty acids: fishy burps and aftertastes. This undesirable effect has even demotivated the otherwise health conscious consumers to avoid taking fish oil supplements altogether.

Until now. GELITA® EC – the first commercially available gelatine product of its kind – allows for the one-step production of enteric capsules that open in the intestine (instead of the stomach, as with traditional gelatine capsules).

! Added Value with GELITA® EC

Most enteric delivery systems are produced by applying an acid-insoluble coating to freshly produced soft capsules. This intensive two-step process adds time and money to the cost of each enteric capsule. Additionally, this coating produces an opaque shell – which is less desired by consumers.

With GELITA® EC, capsules producers now have the ability to create enteric capsules using existing equipment in a one-step process – avoiding additional time and costs – while creating brilliantly clear capsules.

! Proven Enteric Release

GELITA® EC allows manufacturers to produce enteric capsules that adhere to USP (US Pharmacopeia) dissolution parameters. **These parameters include:**

- Capsules will not release more than 10% of the fill during two hours in 37°C simulated gastric fluid
- Capsules will be fully dissolved by 45 min in simulated intestinal fluid

Dissolution results with 20 oblong enteric fish oil capsules produced with GELITA® EC.

- 0% of capsules released the fill in simulated gastric acid (representing the stomach)
- 100% of the capsules showed NO RELEASE (no "leakers")
- 2.15 minutes = average time for capsules to release when pH adjusted to represent intestinal fluids



One Product. One Step.

GELITA® EC Add Value to Consumers and Industry

Benefits of GELITA® EC

- Promotes consumer compliance by eliminating fishy burps often associated with fish oils
- Offers true enteric performance
- Reduces additional time and processing steps associated with historic enteric processes
- Provides an enteric solution with one product and one process
- Although currently available in bovine bone types, GELITA is developing solutions from other raw material sources.

 Patent
Pending



Gelita excludes any warranty and/or liability for any and all claims arising from or in relation with statements made in this brochure and/or the positioning of GELITA® EC, including but not limited to claims based on misleading advertising and/or breach of relevant German and EU legislation. The positioning as well as the packaging, labelling and advertising of the client's products with GELITA® EC falls solely within the responsibility of the client.

GELITA AG · Uferstr. 7 · 69412 Eberbach · Germany · www.gelita.com

GELITA
Improving Quality of Life